




**Klausz Dental  
Laboratories Ltd.**

Working **harder and smarter** for your practice!


www.klauszdentalab.com • klauszdl@uniserve.com • Tel: 416-630-6262 • Fax: 416-630-0163  
40 Kodiak Crescent • Units 11-12 • Toronto, ON • Canada M3J 3G5

## Mandibular Day Orthotics

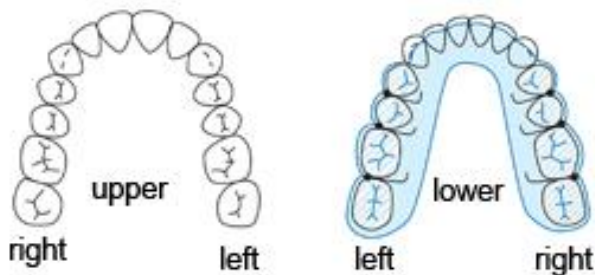


### Mandibular Positioning Appliance

- \* Custom PMT Base or Flexible Thermoplastic Base
- \* Optional ball clasp 3-4 per side
- \* Full coverage over Md posterior teeth
- \* Optional anterior wrap-over incisal coverage
- \* Lingual anatomical ridge 30% of Mx post cusp tips index
- \* No contact of Mx buccal cusps
- \* No protrusive and non-working interferences

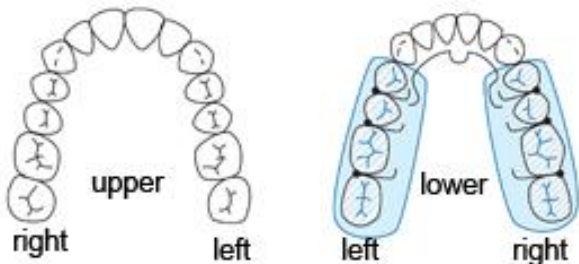


### Compact Mandibular Positioning Appliance

- \* Custom PMT Base
- \* .032 lingual connector wire to replace major connector
- \* Ball clasp 3-4 per side
- \* Full coverage over Md posterior teeth
- \* Lingual anatomical ridge 30% of Mx post cusp tips index
- \* No contact of Mx buccal cusps
- \* No protrusive and non-working interferences

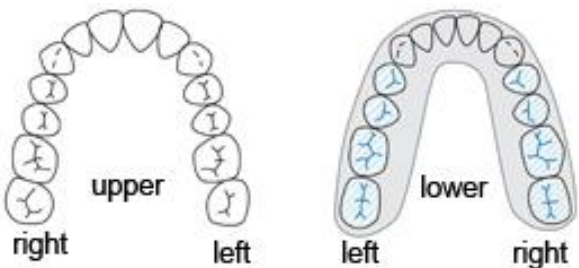

### Mandibular Positioning Appliance (No Ball Claps)

- \* Custom PMT Base
- \* Full coverage over Md posterior teeth
- \* Optional anterior wrap-over incisal coverage
- \* Lingual anatomical ridge 30% of Mx post cusp tips index
- \* No contact of Mx buccal cusps
- \* No protrusive and non-working interferences

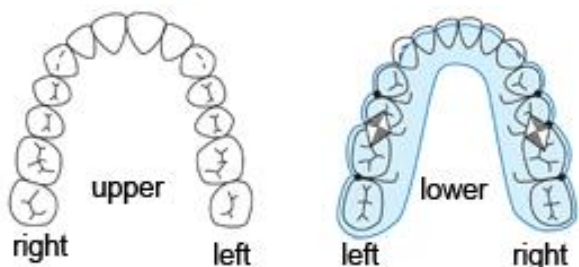



### Pivot Appliance

- \* Custom PMT Base
- \* .032 Ball clasp 3 per side
- \* Full coverage over Md posterior teeth
- \* Optional anterior wrap-over incisal coverage
- \* Pyramid shaped build-up over mesial buccal cusp of lower first molars and occludes against mesiolingual cusp of upper first molars.

***Whatever your requirements, we can custom fabricate an orthotic to meet your specific needs!***





## Maxillary Night Orthotics

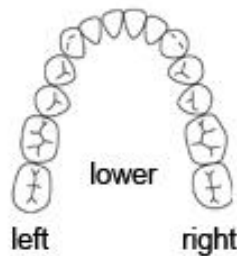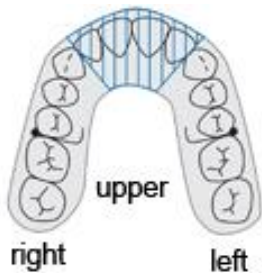



### Anterior Deprogrammer

- \* Custom PMT Base or Flexible Thermoplastic Base
- \* Optional with ball clasp between 5-6's
- \* Anterior acrylic pad contacting lower 4 incisors and extends labially and lingually to patients range of motion
- \* No posterior tooth contact.

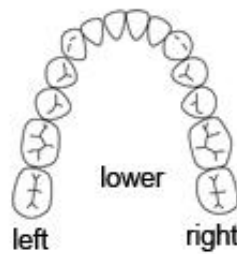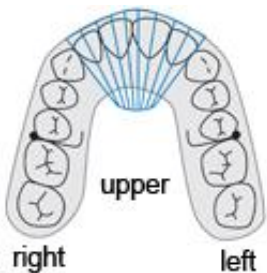



### Anterior Deprogrammer with Lingual Ramp

- \* Custom PMT Base or Flexible Thermoplastic Base
- \* Optional with ball clasp between 5-6's
- \* Anterior acrylic pad contacting lower 4 incisors and extends labially and lingually to patients range of motion
- \* Solid anterior ramp contoured to lingual of lower cuspid to cuspid
- \* No posterior tooth contact.

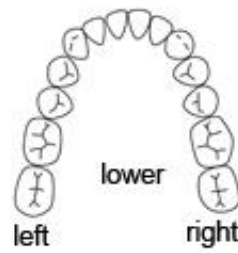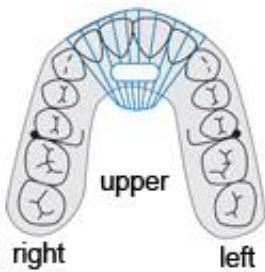

### Anterior Deprogrammer with Open Lingual Ramp

- \* Custom PMT Base or Flexible Thermoplastic Base
- \* Optional with ball clasp between 5-6's
- \* Anterior acrylic pad contacting lower 4 incisors and extends labially and lingually to patients range of motion
- \* Anterior ramp with airway cut out and 2mm frame work contoured to lingual of lower cuspid to cuspid
- \* No posterior tooth contact

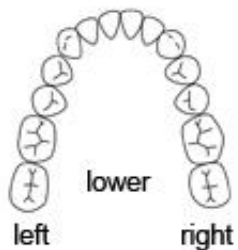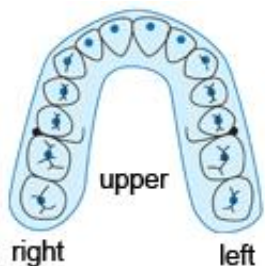

### Decompressor (Flat Plane Splint)

- \* Custom PMT Base or Flexible Thermoplastic Base
- \* Optional with ball clasp between 5-6's
- \* Mx full acrylic coverage and occlusion
- \* Minimum 2mm thickness

*Whatever your requirements, we can custom fabricate an orthotic to meet your specific needs!*